

taylordolman

Guldmann Slings and Accessories



Slings and lifting accessories

Slings and other lifting accessories are a key part of the extensive range of Guldmann equipment for lifting, moving, positioning, training, and mobilising people.

Specially designed slings make it easy, practical and quick for the carer to support the user while being lifted and moved from one point to another, as well as ensuring that the carer can undertake such operations with a minimum of effort and strain.

All Guldmann slings are designed for comfort, dignity and ease of use with both mobile and ceiling hoist systems, and are ideal for moving people to and from beds and wheelchairs, as well as in sitting and lying positions, as required.

Special Guldmann sling designs are available for lifting and moving users with:

- different abilities to participate actively
- different degrees of body control and muscular strength
- different degrees of mobility
- single and double amputations
- special needs

and for both adults and children of many different sizes and weights.

They can also be used as integral parts of training and activity programmes for users. All Guldmann slings are available in a reusable (washable) version and/or in a disposable version to help ensure the best possible hygiene.

Accessories


In addition to slings, Guldmann provides a wide range of specially designed lifting accessories for helping both the user and carer ensure greater mobility with less effort. And they allow the carer to provide the best possible assistance, while working in ergonomically correct positions.




Materials




*Polypropylene
(Disposable)*

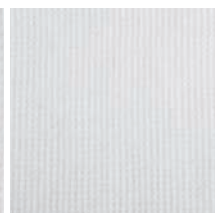


*Polyester/Cotton
(Bariatric)*

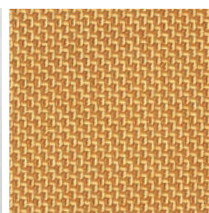

Polyester

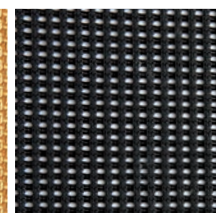

*Net
Polyester*

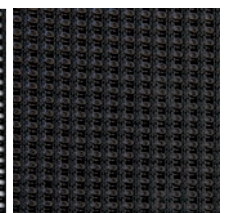

*Spacer
(Comfort)*

ABC Slings

ACTIVE

BASIC

CUSTOM

Slings		Kids	Kids	Kids	XS	S	M	L	XL	XXL	3XL
		4 - 6	6 - 10	10 - 14							
Active Micro	polyester				284031	284041	284051	284061	284071	284081	
Active Micro Plus	polyester	281001	281011	281021	281531	281541	281551	281561	281571	281581	
Active Trainer	polyester		283011	283021	283031	283041	283051	283061	283071	283081	
Basic Basic	polyester	270001	270011	270021	270031	270041	270051	270061	270071	270081	270091 ^{x)}
Basic Basic	net	270102	270112	270122	270133	270143	270153	270163	270173	270183	
Basic Low	polyester				271031	271041	271051	271061	271071	271081	
Basic Low	net				271133	271143	271153	271163	271173	271183	
Basic High	polyester	272001	272011	272021	272031	272041	272051	272061	272071	272081	272091 ^{x)}
Basic High	net	272102	272112	272122	272133	272143	272153	272163	272173	272183	
Basic Shell	net				275032	275042	275052	275062	275072		
Basic Hammock	polyester				274031	274041	274051	274061	274071		
Sit-On	net	295001	295011		295031	295041	295051	295061	295071		
Sit-On High	net	296001	296011		296031	296041	296051	296061	296071		
Sit-On Comfort	spacer	293001	293011		293031	293041	293051	293061	293071		
Sit-On Comfort High	spacer	294001	294011		294031	294041	294051	294061	294071		
Custom Amputee	polyester				290031	290041	290051	290061	290071		
Active Vest Kids ^{xx)}	polyester		283111	283121 ^{xxx)}							
Vest for Standing Shell ^{xx)}	polyester		283511	283521 ^{xxx)}	283531						


Max lifting capacity for all Guldmann ABC slings: 255 kg

- x) Basic Basic and Basic High XXXL can lift up to 500 kg
- xx) Max lifting capacity for Active Vest Kids / Vest for standing Shell: 100 kg
- xxx) Active Vest Kids / Vest for standing Shell: Kids 10-16

Guidelines on choosing the size of a sling

	Kids	Kids	Kids	XS	S	M	L	XL	XXL	3XL
	4 - 6	6 - 10	10 - 14							
H	45 - 60	55 - 70	68 - 80	75 - 82	80 - 87	85 - 92	90 - 97	95 - 102	95 - 102	95 - 102
W	25 - 30	28 - 35	33 - 38	33 - 37	36 - 40	39 - 43	42 - 45	45 - 49	48 - 52	54 - 58

All dimensions in cm

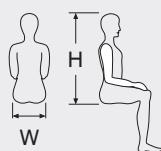

Normal washing at the indicated temperature




Do not use bleaching agent



Tumble-drying at low temperature



Do not iron

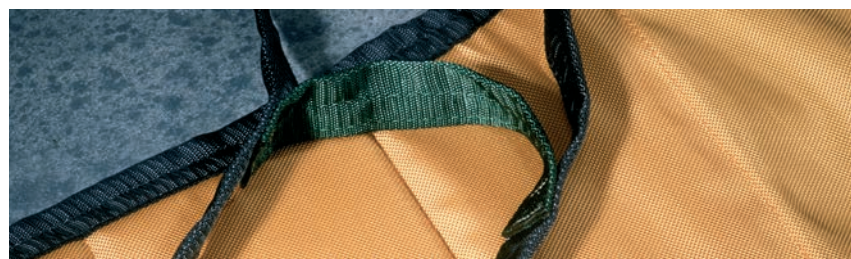

The coloured loops ensure symmetry in the lift and at the same time ease the communication in connection with the lift

ABC Slings



- ① **Coloured loops**
The coloured loops ensure symmetry in the lift
- ② **Material**
Dirt repellent polyester characterized by its strength and elasticity
- ③ **Auxiliary straps**
It is possible to use auxiliary straps

- ④ **Padding on legs**
All slings with divided leg supports have inset foam pads on legs. Easier to put on – smoother on skin
- ⑤ **Centre of sling marking**
- ⑥ **Information-label**
Indicates all that the carer needs to know

- ⑦ **Label with size indication**
Makes it easy and simple to find the correct size of sling
- ⑧ **Padding in back section**
For added comfort on models with small back surface
- ⑨ **Positioning pocket**
All Basic slings are provided with a pocket to ease correct positioning of sling

Fit
Improved fit and more sizes to choose from for added comfort and safety

Strength
Max. capacity increased to 255/500 kg

Washing
Easy handling, the sling can be washed at 80°C and tumble-dried at max. 50°C



ABC Slings – Active



1 Active Trainer

Function

Provides support around the chest and shoulder areas. Detachable leg straps give support around the hip area and prevent the user from sliding out of the sling.

Suitable for

- persons with reduced body balance, but must be able to bear weight on their legs
- walking training • supported standing and upright postures • early mobilization.

..... 2830xx

2 Active Micro

Function

Supports the upper part of the body (from the pelvis to under the shoulders) around the chest and upper thigh. The chest straps tighten around the chest.

Suitable for

- lifting persons with head control, but with a reduced muscular strength of upper part of the body and above hip joint and thighbone
- for some persons with extension spasms
- in situations where the sling should be put on and taken off easily • lifting to and from a seated position • situations where a large opening at the bottom is required, such as access for toileting and hygiene.

..... 2840xx

3 Active Micro Plus

Function

Supports the upper part of the body (from the pelvis to under the shoulder) around the chest and upper thigh. Support belt can be adjusted.

Suitable for

- lifting persons with head control, but with a reduced muscular strength of upper part of the body and above hip joint and thighbone
- in situations where the sling should be put on and taken off easily • lifting to and from a seated position • situations where a large opening at the bottom is required, such as access for toileting and hygiene.

..... 2815xx



ABC Slings – Basic

4



5



6



4

Basic Basic

Function

Supports the whole body excluding the head.

Suitable for

- lifting persons with head control and with or without body control
- lifting to and from a seated or half-seated position (bed and wheelchair etc.)
- lifting and transfers from all patient furniture surfaces.

..... 270xxx

5

Basic Low

Function

Supports the body up to shoulder height.

Suitable for

- lifting persons with head and body control
- lifting to and from a seated position, (chair and bed etc.)
- lifting and transfers from all patient furniture surfaces.

..... 271xxx

6

Basic High

Function

Supports the whole body including the head.

Suitable for

- lifting persons with little control over their head and body
- lifting to and from a lying position (bed and floor etc.)
- lifting to and from a seated position
- lifting and transfers from all patient furniture surfaces and the floor.

..... 272xxx



ABC Slings – Basic



7
Basic Shell

Function
Supports the whole body and the head, but allows the arms to be outside the sling.

Suitable for
• lifting persons with little control over their head and body • lifting to and from a lying position (bed and floor etc.) • lifting to and from a seated position – letting the sling stay in the wheelchair.
..... 2750xx



8
Basic Hammock

Function
Supports the whole body apart from the head. If head support is required, the head support can be used in conjunction with the sling. Adds extra support around the hip and thigh area. Can be used as divided legs or as a cross over position, this also provides good support for the amputee (amputee users **must** be assessed for this sling).

Suitable for
• lifting persons with little control over their body • lifting to and from a lying position (bed and floor etc.) • lifting to and from a seated position • persons who need added support around the hip and thighs.
..... 2740xx

ABC Slings – Custom

9



10



11



9

Sit-On

Function

Supports the thighbone and back. The user sits on the sling.

Suitable for

Persons who remain on the sling for long periods due to high weight or customized wheelchair.

Other features

- helps reduce the risk of sliding out of wheelchair
- material prevents heat generation between body and sling
- sling is fitted on the person while lying, but is used when lifting to and from a seated position in bed or wheelchair.

..... 2950xx

10

Sit-On High

Function

Supports the thighbone, back and head. The user sits on the sling.

Suitable for

Persons who remain on the sling for long periods due to high weight or customized wheelchair.

Other features

- helps reduce the risk of sliding out of wheelchair
- material prevents heat generation between body and sling
- sling is fitted on the person while lying, and is used when lifting to and from a lying position.

..... 2960xx

11

Custom Amputee

Function

Supports the whole body up to shoulder height and the pelvis region.

Suitable for

- lifting both single and double (above knee joint) amputees
- lifting double and high amputees
- lifting from a lying position. Head support should be used to support head.

..... 2900xx



ABC Slings – Custom

12



13



12

Sit-On Comfort

Function

Supports the thighbone and back. The user sits on the sling. The design of the sling prevents the thighbone from an invert rotation when hoisting. Added comfort due to the Spacer material.

Suitable for

Persons who remain on the sling for long periods due to high weight or customized wheelchair.

Other features

- helps reduce the risk of sliding out of wheelchair
- material prevents heat generation between body and sling
- sling is fitted on the person while lying, but is used when lifting to and from a seated position in bed or wheelchair.

..... 2930xx

13

Sit-On Comfort High

Function

Supports the thighbone, back and head. The adjustable head support can be folded down when user is left in sitting position on the sling. The design of the sling prevents the thighbone from an invert rotation when hoisting. Added comfort due to the Spacer material.

Suitable for

Persons who remain on the sling for long periods due to high weight or customized wheelchair.

Other features

- helps reduce the risk of sliding out of wheelchair
- material prevents heat generation between body and sling
- sling is fitted on the person while lying, and is used when lifting to and from a lying position.

..... 2940xx

ABC Slings – Kids

14



15



16



14

Basic Basic – Kids

Function

Supports the whole body excluding the head. The inside belt offers an extra support of the body.

Suitable for

- lifting kids with head control and with or without body control
- lifting to and from a seated or half-seated position (bed and wheelchair etc.).

..... 270xxx

15

Basic High – Kids

Function

Supports the whole body including the head. The inside belt offers an extra support of the body.

Suitable for

- lifting kids with little control over their head and body
- lifting to and from a lying position (bed, changing bench and floor etc.)
- lifting to and from a seated position.

..... 272xxx

16

Active Vest – Kids

Function

Supports the body and pelvis. The width can be adjusted.

Suitable for

- kids with reduced body balance
- walking training
- supports standing and upright postures
- playing.

..... 2831xx

ABC Slings – Kids

17



18



19



17

Sit-On Comfort – Kids

Function

Supports the thighbone and back. The user sits on the sling. The design of the sling prevents the thighbone from an invert rotation when hoisting. Added comfort due to the Spacer material.

Suitable for

Kids who remain on the sling for long periods due to customized wheelchair, or need hoisting regularly.

Other features

- helps reduce the risk of sliding out of wheelchair
- material prevents heat generation between body and sling
- sling is fitted on the child while lying, but is used when lifting to and from a seated position in bed or wheelchair

..... 2930xx

18

Sit-On Comfort High – Kids

Function

Supports the thighbone, back and head. The adjustable head support can be folded down when user is left in sitting position on the sling. The design of the sling prevents the thighbone from an invert rotation when hoisting. Added comfort due to the Spacer material.

Suitable for

Kids who remain on the sling for long periods due to high weight or customized wheelchair.

Other features

- helps reduce the risk of sliding out of wheelchair
- material prevents heat generation between body and sling
- sling is fitted on the child while lying, and is used when lifting to and from a lying position.

..... 2940xx

19

Vest for Standing Shell – Kids

Function

Supports the head and back

Suitable for

- lifting the child from a lying to a standing position
- used in conjunction with a standing shell to assist into a standing frame or table.

..... 2835xx



Slings – Bariatric

20



20

Basic High Bariatric, 500 kg

Function

Supports the whole body including the head.

Suitable for

- lifting persons with little control over their head and body
- lifting to and from a lying position (bed and floor etc.)
- lifting to and from a seated position
- lifting and transfers from all patient furniture surfaces and the floor.

4XL standard 2720914

Slings – Bariatric



21



22



23



21

Twin Turner Bariatric
 Makes it possible to roll a patient on both sides without removing the sling. Enables placing of sling under patients without having to roll or move them manually. Smaller surface area of material for this sling allows access to patient's skin for inspection, treatment and care.
 28760

22

Leg sling
 Used when legs need lifting and support. Possible activities are dressing, apply therapeutic sock/hose, access to hygiene, support of leg to apply ROM devices.
Packed in boxes of 10 slings.
 28650

23

Repositioning sling Bariatric, 500 kg
 A sling for use as part of the linen make-up of the patient bed; leave positioned underneath patient, ready to use for care tasks. Suitable for: rolling to sidelying, repositioning patient up in bed, sitting forward head and trunk, supine lateral transfers to another surface such as beds + stretchers, and full lift in supine off of bed to zero bed scale, stripe linen, or other.
Packed in boxes of 5 slings.
 284656

Slings – Bariatric

24



24

Gait Trainer Bariatric, 500 kg

Function

Supports the pelvis area and pannus.

Suitable for

bariatric patients. Will support the patient in a standing position, maximising comfort and support.

.....283100





Disposable slings

25



25

Disposable High II, 205 kg

Function

Supports the whole body including the head. This is a single patient use sling that tolerates water – dispose of after discharge. Beneficial for infection control program. Same design as Guldmann Basic High sling.

Suitable for

- lifting persons without head and body control
- lifting to and from a lying position (bed, floor, etc)
- lifting to and from a seated position.

Packed in boxes of 10 slings.

Material

Polypropylene

XS/S	283641
S/M	283651
M/L	283661
L/XL	283671
XL/XXL	283681

26

Disposable High Bariatric, 350 kg

Function

Supports the whole body including the head. This is a single patient use sling that tolerates water – dispose of after discharge. Beneficial for infection control program. Same design as Guldmann Basic High sling.

Suitable for

- lifting persons without head and body control
- lifting to and from a lying position (bed, floor, etc)
- lifting to and from a seated position.

Packed in boxes of 10 slings.

Material

Polyester/cotton

L/XL	283672
XL/XXL	283682
XXL/XXXL	283692

27

Disposable High II – Kids, 205 kg

Function

Supports the whole body including the head. This is a single patient use sling that tolerates water – dispose of after discharge. Beneficial for infection control program. Same design as Guldmann Basic High sling.

Suitable for

- lifting kids without head and body control
- lifting to and from a lying position (bed, floor, etc)
- lifting to and from a seated position.

Packed in boxes of 10 slings.

Material

Polypropylene

Kids 6-10	283611
-----------	--------



Disposable slings

28



28

**Disposable
Repositioning sling, 500 kg**

Function

A sling for use as part of the linen make-up of the patient bed; leave positioned underneath patient, ready to use for care tasks. This is a single patient use sling that tolerates water – dispose of after discharge. Beneficial for infection control program. Same design as Guldmann Repositioning sling.

Suitable for

Rolling to sidelying, repositioning patient up in bed, sitting forward head and trunk, supine lateral transfers to another surface such as beds + stretchers, and full lift in supine off of bed to zero bed scale, change linen, or other.
Packed in boxes of 10 slings.

Material

Polyester
..... 284655

29



29

**Disposable
Horizontal sling standard, 350 kg**

Function

Supports the whole body.

Suitable for

Lifting a patient in a horizontal position.
Packed in boxes of 10 slings.

Material

Polyester
..... 284631

Disposable slings

30
Disposable
Twin Turner II

Function
 To be used when turning a user into lateral position. This is a single patient use sling that tolerates water. Enables placing of sling under the patient without having to roll or move them. Smaller surface area of material for this sling allows access to patient's skin for inspection, treatment and care. Same design as Guldmann Turner.
Packed in boxes of 10 slings.

Material
 Polypropylene
 Regular 287501
 Large 287511



31
Disposable
Leg sling II

Function
 Used when legs need lifting and support. Possible activities are dressing, apply therapeutic sock/hose, access to hygiene, support of leg to apply ROM devices. Dispose of after discharge. Beneficial for infection control program. Same design as Guldmann Leg sling.
Packed in boxes of 10 slings.

Material
 Polypropylene
 286501



32
Disposable
Multi Support sling

Function
 Supports lifting part of the body, upper body, pelvis, legs or arms.

Suitable for
 Situations where parts of the body need to be lifted or supported – e.g. for placing X-RAY console.
 284223





Turning and handling options



33
Positioning sling

The positioning sling system makes it possible to easily turn the patient to one side or to a prone position on e.g. the OR table. Using the sling together with the hoist prevents shear force and friction on the patient's skin and improves the working environment for the care staff.

Side positioning sling 284226

Prone positioning package:

2 x slings

2 x Multi Support slings

1 x security strap

. 284225

Turning and handling options



35
Leg sling
 Used when legs need lifting and support. Possible activities are dressing, apply therapeutic sock/hose, access to hygiene, support of leg to apply ROM devices.
Packed in boxes of 10 slings.
 28650



36
Repositioning sling, 500 kg
 A sling for use as part of the linen make-up of the patient bed; leave positioned underneath patient, ready to use for care tasks. Suitable for: rolling to sidelying, repositioning patient up in bed, sitting forward head and trunk, supine lateral transfers to another surface such as beds + stretchers, and full lift in supine off of bed to zero bed scale, stripe linen, or other.
Packed in boxes of 10 slings.
 Standard 28465
 Short (photo) 284653
 Net. 284650

Turning and handling options



37

Turner

Enables placing of sling under patient without having to roll or move them. Smaller surface area of material for this sling allows access to patient's skin for inspection, treatment and care 28700

38

Twin Turner

Makes it possible to roll a patient on both sides without removing the sling. Enables placing of sling under patients without having to roll or move them. Smaller surface area of material for this sling allows access to patient's skin for inspection, treatment and care.

Regular 28750
Large 28751

39

Multi Support sling

Function

Supports lifting part of the body, upper body, pelvis, legs or arms.

Suitable for

Situations where parts of the body need to be lifted or supported – e.g. for placing X-RAY console. 28467



Stretcher devices

40



41



40

Horizontal lifter, foldable, 250 kg 2 strap version

Suitable when a horizontal and stable lift is required, with weight adjustment. To be used with GH3 Twin. Can be used with either a complete lifting sheet or straps.

- Horizontal lifter
- incl. one lifting sheet 28444
- Lifting sheet 284485
- Set of lifting straps incl. stiffener 28448

41

Horizontal lifter, foldable, 250 kg

Suitable when a horizontal and stable lift is required, with weight adjustment. To be used with GH3. Can be used with either a complete lifting sheet or straps.

- Horizontal lifter
- incl. one lifting sheet 28445
- Lifting sheet 284485
- Set of lifting straps incl. stiffener 28448

42

Trolley for horizontal lifter

. 28443



Stretcher devices



43



46



47



44



45



48

43
Horizontal lifting support, 350 kg
 Horizontal lifting support with infinitely variable weight adjustment. Can be used with GH3 ceiling hoists or with Guldmann mobile lifter.
 28456
 Horizontal sling, standard 28463

44
Horizontal lifting support, 350 kg
 Horizontal lifting support with 5 different lifting positions. Can be used with GH3 ceiling hoists or with Guldmann mobile lifter.
 28466
 Horizontal sling, standard 28463

45
Mortuary lifting straps
 Ease lifting in mortuary. Can take larger corpses / bodies. The Individual straps allow easy access.
 Standard incl. stiffener 284640

46
Bath chair – rustproof, 150 kg
 Suitable when a stable lift into a swimming pool is required.
 14600

47
Bath stretcher – rustproof, 150 kg
 Suitable when a stable and horizontal lift into a swimming pool is required.
 14610

48
Horizontal lifter BS, 150 kg
 A horizontal lifter especially developed and designed to overcome the many moving and handling challenges that mortuary staff faces every day.
 28425

Lifting accessories



49
Digital Scale for Ceiling hoist or Mobile lifter
 Battery operated digital scale.

(A) Max lifting capacity 350 kg
 Approved for medical use. 11346

(B) Max lifting capacity 250 kg 113461


51
Head support
 Supports the head without stimulating the neck reflexes. Suitable for persons with no head control. 28325

52
Extension straps
 Used in situations to achieve various sitting positions or where different lifting heights are required. *Set of two pcs.*

20 cm 28470
 30 cm 28472

50
Sheepskin insets
 Suitable for persons who need more gentle support of the thighbone. 28315

Lifting hanger options

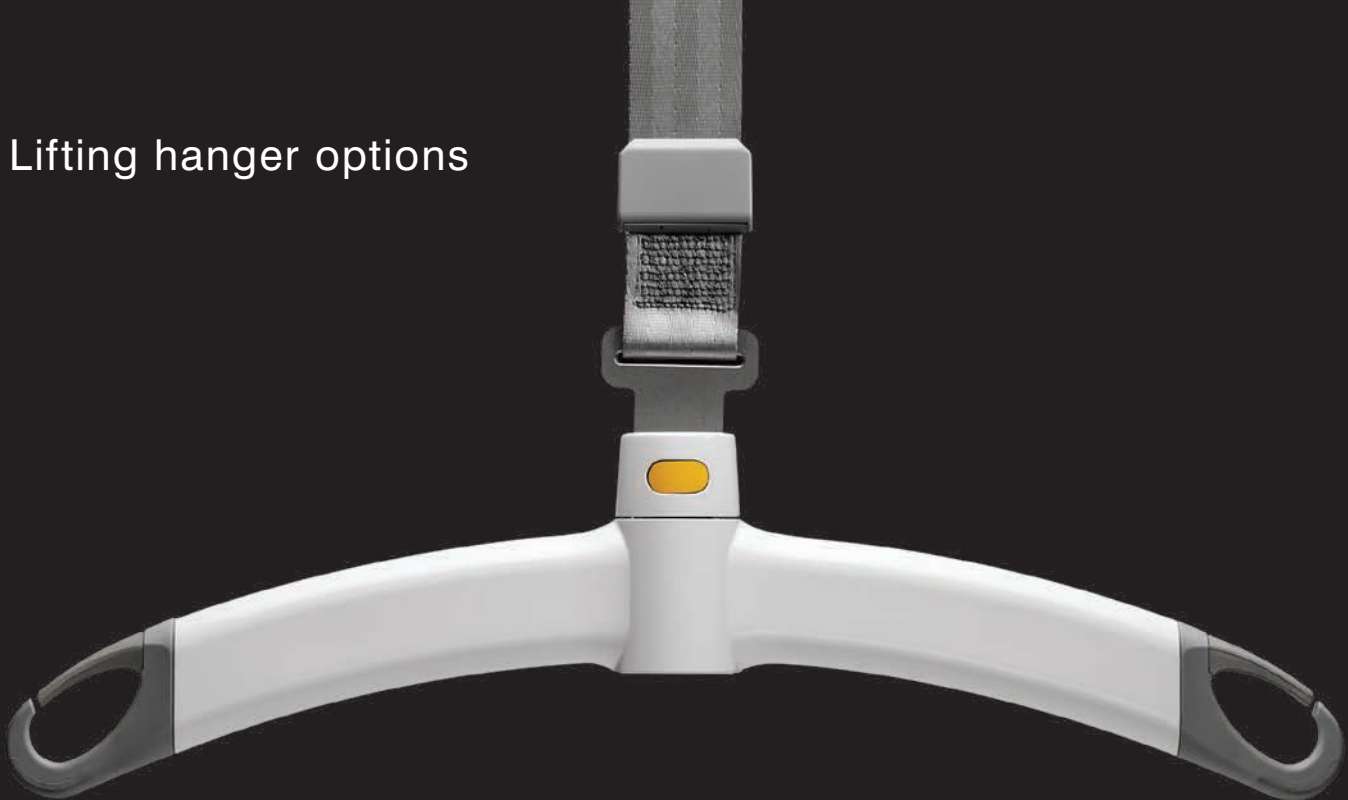

53



53
GH Lifting hanger, medium, 350 kg
 Suitable for persons using a lifting sling in sizes M-XXXL.

Function
 To be used in connection with uncomplicated lifts and when using standard slings.

..... 550316

54



54
GH Lifting hanger, small, 350 kg
 Suitable for persons using a lifting sling in sizes XS and S.

Function
 To be used in connection with uncomplicated lifts and when using standard slings. The lifting hanger is recommended when using lifting slings size small and down. Can also be used with the Active Micro Plus lifting sling.

..... 550057

55




55
GH Lifting hanger, x-small, 350 kg
 Suitable for persons using a lifting sling in XS and Kid Sizes.

Function
 To be used in connection with uncomplicated lifts and when using standard slings. The lifting hanger is recommended when using the "kids" lifting slings.

..... 550058

Lifting hanger options

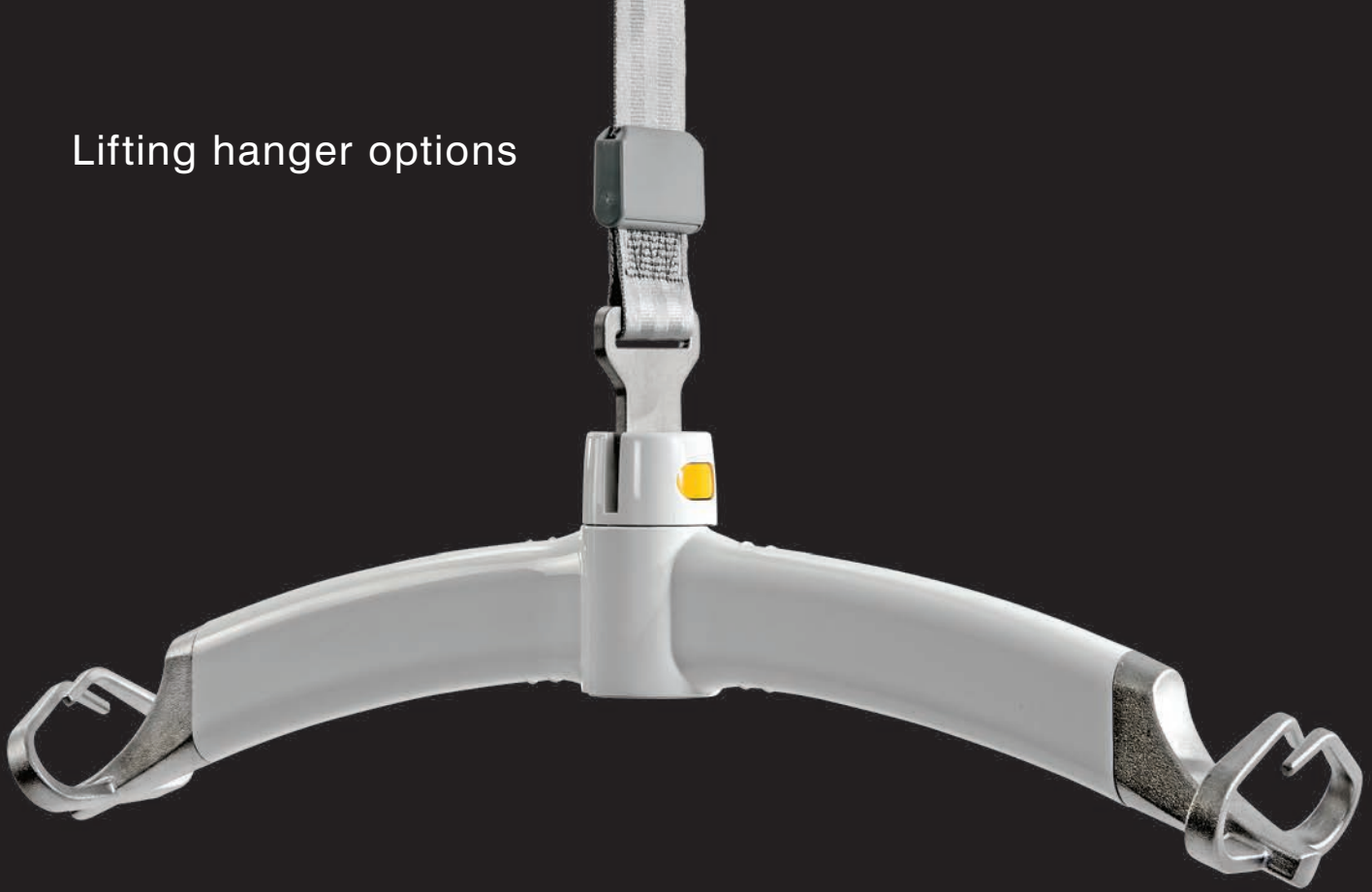

56



56
GH lifting hanger, 4 attachment points medium, 350 kg

Suitable for persons using a lifting sling in sizes M-XXXL.

Function

To be used when you need a small distance between the 4 lifting points and when you are working with slings where you need more lifting straps in the same lifting hook.

..... 554258

57



57
GH lifting hanger, 4 attachment points small, 350 kg

Suitable for persons using a lifting sling in sizes XS and S.

Function

To be used when you need a small distance between the 4 lifting points and when you are working with slings where you need more lifting straps in the same lifting hook.

..... 554259

58



58
GH lifting hanger, 4 attachment points x-small, 350 kg

Suitable for persons using a lifting sling in XS and Kid Sizes.

Function

To be used when you need a small distance between the 4 lifting points and when you are working with slings where you need more lifting straps in the same lifting hook.

..... 554260

Lifting hanger options



59



60



61

59

GH Cross hanger, 500 kg

Suitable when larger spacing between back and leg straps of the lifting sling is required. To be used when lifting with GH3+ and GH3 Twin ceiling hoist.

..... 550800

60

GH Cross hanger, 250 kg

Suitable when larger spacing between back and leg straps of the lifting sling is required. To be used when lifting with GH3+ and GH3 Twin ceiling hoist.

..... 554904

61

H Hanger, 350 kg, GH3

Suitable when lifting persons where a large spacing between back and leg straps of the lifting sling is required, e.g. when lifting persons with a large abdomen/ stomach. The H hanger is available with side hangers in 3 different sizes.

H Hanger – Centre hanger 550443

(Consists of 1 hanger 550316 and 2 Swing kit 554903)

Note: In addition you need to order 2 hangers of either x-small (550058), small (550057) or standard (550316)

Lifting hanger options



62
Oval Hook
 21160

63
Swing kit
 554903

64
Parking space for GH lifting hanger
 550531
 Base for parking space for
 GH lifting hanger, 4 attachment points
 551329

65
Connecting bar
 Optional for GH3 Twin so diverse
 hangers can be utilized
 550544

taylordolman
 tel: 0845 838 0600
 email: info@taylordolman.com
www.taylordolman.com