

Etac Molift Quick Raiser 1

Compact design and excellent maneuverability even in narrow spaces.



Compact and easy

By virtue of its compact design and low weight, Etac Molift Quick Raiser 1 is stable and easy to maneuver. The unique wheel position and near perfect weight distribution of the base provides a small turning radius with easy maneuverability even in narrow spaces.

Natural pattern of movement

The inclined straight column imitates the natural pattern of movement and fulfills the lifting needs of users with some degree of weight bearing capacity. Molift Quick Raiser 1 has a safe working load of 160 kg (352 lbs.)

The extra low base, 7 cm (2.8") and narrow front, 34 cm (13.5") fits well around toilets, under most beds, chairs and wheelchairs.

Find detailed information at page 12–13.

Features



Extra low base
Extra low base with narrow front for easy access.



Integrated service software
Automatic service notification system determines when the lifter requires service and maintenance.



Compact and easy to maneuver
Etac Molift Quick Raiser 1 has a low weight and is stable and easy to maneuver.



Multiple adjustable knee pads
The knee pads are height adjustable and can be turned around for better fitting and to give alternative support for the legs.

taylordolman



Etac Molift Quick Raiser 2

Sit-to-stand passive lift with electrically spreadable leg base.



Stable and easy to maneuver

Etac Molift Quick Raiser 2 with electrically spreadable leg base is compact, stable and easy to maneuver. The unique wheel position and near perfect weight distribution of the base provides a good turning radius and easy maneuverability even in narrow spaces.

Fits almost everywhere

The inclined straight column imitates the natural pattern of movement and fulfills the lifting needs of users with some degree of weight bearing capacity. Molift Quick Raiser 2 has a safe working load of 160 kg (352 lbs)

The low base, 10 cm (4") fits under most beds and chairs whilst the electrically spreadable legs enables accessibility almost everywhere.

Find detailed information at page 12-13.

Features



Electrically spreadable legs

The perfect alternative when the lifter is used together with chairs with no accessibility underneath, such as electric wheelchairs and recliners.



Natural pattern of movement

Inclined straight column imitates a natural pattern of movement.



Multiple adjustable knee pads

The knee pads are height adjustable & can be turned around for better fitting & to give alternative support for the legs.



Integrated service software

Automatic service notification system determines when the lifter requires service and maintenance.

Etac Molift Quick Raiser 2+

When a higher safe working load is needed.



High safe working load

Etac Molift Quick Raiser 2+ has been adapted for higher safe working load, up to 200 kg (440 lbs). It has a reinforced column and suspension and electrically spreadable leg base. The unique wheel position and perfect weight distribution of the base provides a good turning radius and easy maneuverability even in narrow spaces.

Natural pattern of movement

The inclined straight column imitates the natural pattern of movement and fulfills the lifting needs of users with some degree of weight bearing capacity.

The low base, 10 cm (4") fits under most beds and chairs whilst the electrically spreadable legs enables accessibility almost everywhere.

Find detailed information at page 12–13.

Features



Increased safe working load

Etac Molift Quick Raiser 2+ has a safe working load of 200 kg (440 lbs)



Multiple adjustable knee pads

The knee pads are height adjustable and can be turned around for better fitting and to give alternative support for the legs.



Integrated service software

Automatic service notification system determines when the lifter requires service and maintenance.



Natural pattern of movement

Inclined straight column imitates a natural pattern of movement.

Features – Molift Quick Raiser



Environmentally friendly NiMH battery

All Molift lifts are equipped with an environmentally friendly NiMH-battery.



Battery indicator

All Molift lifts have an integrated battery indicator that provides audio and visual notification.



Unique wheel positioning

The lightweight construction combined with the unique wheel positioning, provides exceptional maneuverability.



Safety features

Numerous safety features are standard, including emergency stop, overload protection, electrical and mechanical emergency lowering.



Power-save function

The electronics has a power-save function which makes the system go into sleep-mode after a period with no activity.

Accessories



Molift Quick Raiser adjustable footboard

Footboard for adjusting the height of Molift Quick Raiser foot base, ideal for users of smaller height or children. Has a SWL of 75kg (165 lbs)
Item no: 0451000



Molift Quick Raiser Sandals


Designed to ensure the feet are held in place during lifting.
Sandals Small: 0452100
Sandals Medium: 0452200
Sandals Large: 0452300



Extra soft pads for knee support

For patients with extra sensitive legs
Soft pads Small: 0455111
Soft pads Large: 0455114

Options

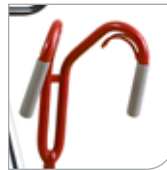

Molift Quick Raiser Arm V-Shape

Is the original lifting arm. It is easy to use and fits most users.
Item no 0458107



Molift Quick Raiser Arm U-shape

The newly developed lifting arm has more space between the attachment points for the sling.
Item no 0458010



Molift Quick Raiser Knee Support Small


Item no 0455000



Molift Quick Raiser Knee Support Large

Item no 0455001

Dimensions

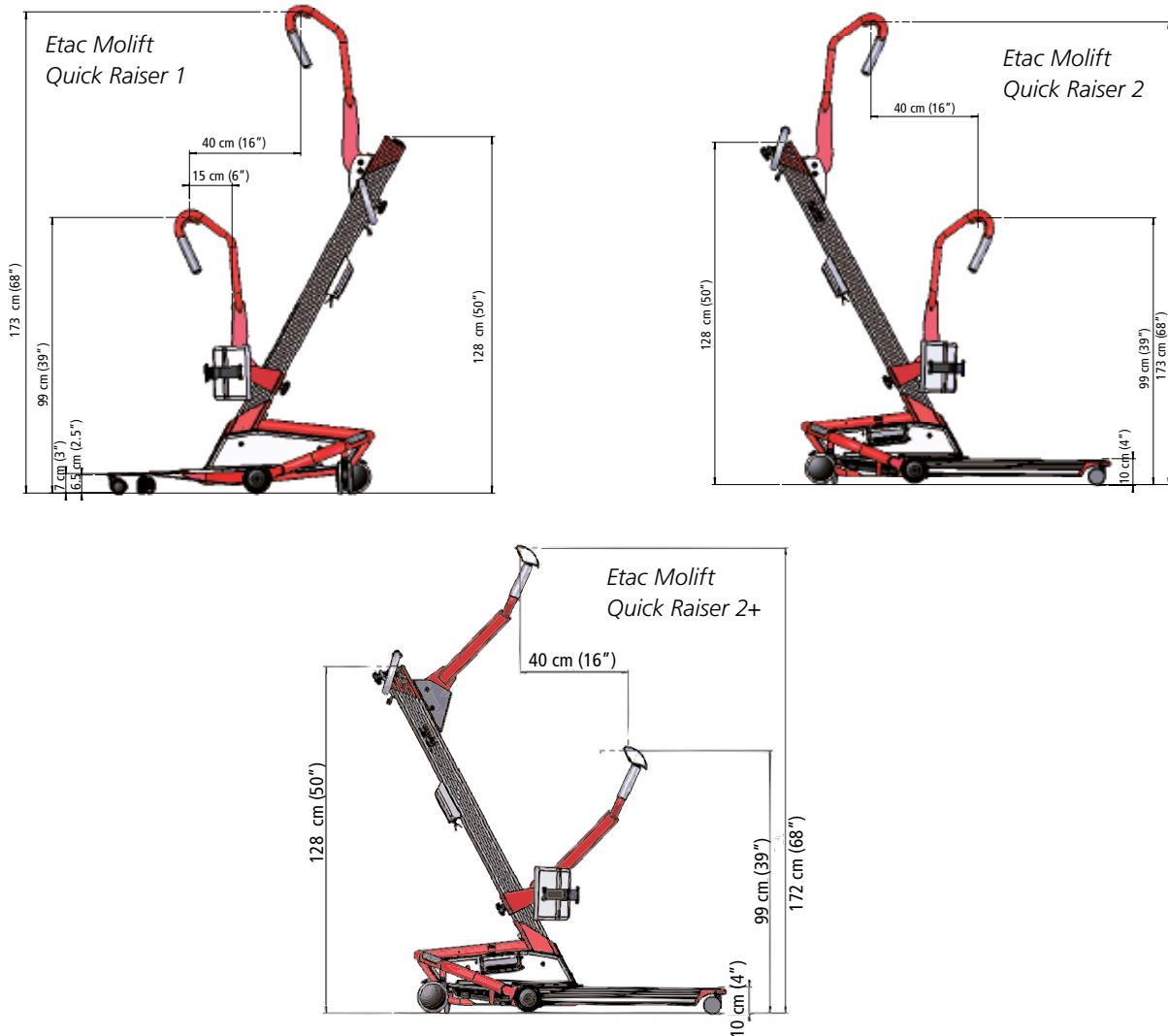

Technical facts Etac Molift Quick Raiser

Description

Molift QR1, V-shape Arm	04100
Molift QR1, U-shape Arm	04150
Molift QR2, V-shape Arm	23160
Molift QR2, U-shape Arm	04200
Molift QR2+, U-shape Arm	04400

Service Software

Included

Material

Aluminium, steel

Total Weight

QR1:	29 kg (64 lbs)
QR2:	36 kg (79.5 lbs)
QR2+:	37 kg (81.5 lbs)

Length

QR1:	95.5 cm (37.6")
QR2 and QR2+:	115 cm (47.2")

Leg height

QR1:	7 cm (2.8")
QR2:	10 cm (3.9")
QR2+:	10 cm (3.9")

Chassis width

61 cm (24.0")

Battery

NiMH 14.4 V - 2.2 Ah

Battery charge time

ca. 3 hours

Lifting speed

5 cm / sec	at 75 kg
(2" / sec)	at 165 lbs

Rated performance

40 lifts at: 75 kg and 50 cm / (165 lbs and 19.5")

Standard item model includes Chassis, Column, Knee Support and respective V- or U-shape Arm

QR1 and 2 QR 2+



Cleaning instructions

Use pH-neutral cleaning agents. If disinfection is required, use isopropyl alcohol. Do not use abrasive or corrosive substances for cleaning.

