



# ROPOX SHOWER/CHANGING BED

### SAFE AND EASY TO USE

The shower/changing bed gives the user optimum comfort and safety, and the helper gets the best possible working position. This provides security and satisfaction for both parties.

The bed is made of double-sided slats. One side of the slat curves slightly downwards ensuring a safe position for the person bathing. The other side curves slightly upwards making it easy for a helper to turn the user on the side when using the bed for changing clothes.

It is easy to adjust the height of the bed with the hand control, and the bed can quickly be adjusted to the best position, e.g. for a wheelchair so transferring a user onto the bed is smooth and easy.

The shower/changing bed is cleaning-friendly, and can be folded up when not in use to save space in the bathroom.








# ROPOX SHOWER/CHANGING BED


Height adjustment:


nt: 70 cm (30-100 cm) above floor

Maximum load: 220 kg cf. DS/EN 17966:2016

#### VARIANT/MODEL

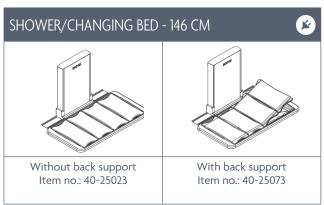

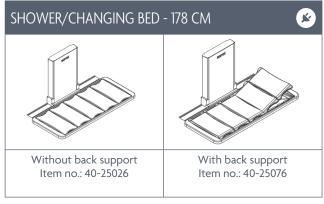

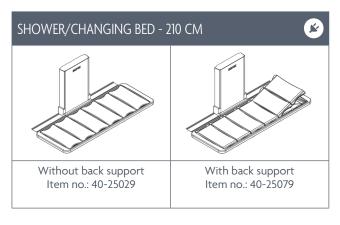



The shower/changing bed includes an electrical hand control for easy height adjustment.

The lying area is 18 cm shorter than the length of the shower/changing bed.

Bed and end guards can be mounted if needed. The transparent bed guard keeps the water spillage to a minimum.

#### **ACCESSORIES**

Back support Item no.: 40-25031



End guard Item no.: 40-25038



Bed guard Item no.: 40-25036



Headrest Item no.: 40-25090



NB! Find more information about measurements, technical data, installation interfaces and much more under the product and in the instruction manuals on our website.